PBM Ortho

## PBM Ortho USER MANUAL

### **INSIDE YOUR PACKAGE:**

PBM Healing Ortho:

- 1 detachable mouthpiece
- 1 battery pack
- 1 power cable connecting battery to mouthpiece
- 1 charging cable for battery
- 1 power adapter
- User manual

---


Before your first use, ensure the battery pack is fully charged.

## HOW TO CHARGE THE BATTERY PACK

- To charge, insert the charging cable into the USB port of the battery pack and connect to the power adapter.
- When charging the 4 white lights near to the charging cable connection will blink.
- Full charge is indicated when the white lights become solid. A full charge takes up to 6 hours.
- A fully charged device is ready to be used for approximately eight sessions.
- As the device is powered by a lithium battery, the battery capacity will gradually decrease with long term usage.

## HOW TO USE THE MOUTHPIECE

1. Connect the mouthpiece to the battery pack via the power cable.
2. Place the mouthpiece in your upper jaw inside your mouth, make sure that it is centered in your jaw (use front teeth as reference).
3. Gently hold the mouthpiece in place with light biting pressure.
4. Press the on-off switch located at the middle of the power cable, and you will see a blue light go on and the treatment starts.
5. After 4 minutes you will hear a beep sound and the blue light will turn off. Treatment for the upper teeth is complete.
6. Flip the mouthpiece, turn the device on again and repeat the above steps on the lower jaw.

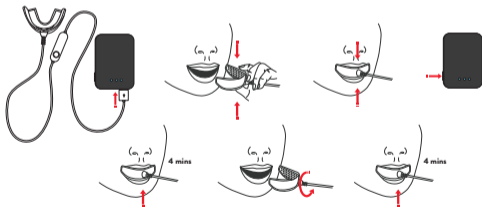TOTAL TREATMENT TIME:

Total - **8** minutes - **4** minutes upper arch and **4** minutes lower arch

A warm and relaxing feeling might be felt in your jaw/gingival area after a while of usage, this is a natural result of light therapy.

## PAUSING USAGE

- The usage can be paused simply by removing the mouthpiece from your mouth.
- You can resume use by putting the mouthpiece back into your mouth.
- If you have paused the device for more than 30 seconds, it would be considered an unsuccessful session and the whole session for that specific jaw would need to be repeated.
- To restart the session, turn off the power and turn it back on. Reinsert the mouthpiece back into the mouth, press the “on-off switch” again, and repeat the 4 minutes for the arch previously being treated.

## COMPLIANCE RECOMMENDATIONS

It is recommended to select the same time every day to do your treatment. It actually takes about 3 weeks to make a habit, so patiently keep using it regularly.

### **Suggestions on making it a habit:**

- Set up PBM healing ortho next to your bed so you can use it the first thing upon rising or just before sleeping.
  - Or place it by your computer or tv and use the PBM healing ortho when using the computer or watching tv.
- the PBM Healing ortho can be used with or without the aligners.

## **MAINTENANCE**

### **Cleaning:**

- Detach the mouthpiece from the battery pack.
- Wash it gently, using clean warm water or the aligner cleaner recommended by PBM Healing.
- To dry, lightly rub the mouthpiece with a clean dry towel and let it air dry. Do not shake or pat down hard on the mouthpiece to dry it.

### **Storage & Handling:**

- Store all components separately. Unplug the battery pack from the charging cable and detach the mouthpiece from the battery pack.
- The PBM Ortho device should be stored in a cool, dry place away from direct sunlight. Avoid storing it in locations where it might be exposed to elevated temperatures.
- It is recommended to store in the original packaging case when it is not intended to be used for a while.
- Special care should be taken, particularly when traveling with the device. This will prevent probable damage and avoid battery depletion.
- Rough handling may cause damage to the device. Stop usage if damage is suspected.

### **Shelf Life:**

- The PBM Ortho device should last for up to two years if used daily with care and within the ranges recommended by the manufacturer.
- The PBM Ortho device includes a lithium battery that will lose charge over time if not recharged. The device should be fully charged prior to the first usage.
- It is highly recommended that the first usage be conducted within two months from product delivery.

### **Part Replacement:**

- All the components are not replaceable or serviceable (including lithium battery).
- Normally, during the usage, components should not need any replacements.
- Wear marks, especially bite marks that become present in the mouthpieces over time of usage, are normal and there is no need to replace the device.
- If the device tears, or the internal surface of the mouthpiece becomes exposed, stop use immediately and contact your doctor or PBM Healing.

## **Benefits of PHOTOBIMODULATION**

PBM Ortho shortens the whole orthodontic treatment process in half and reduces pain associated with orthodontic treatment.

It is clinically and scientifically proven that intense near-infrared light provides a therapeutic benefit. Infrared light energy emitted from the PBM Ortho device into the gingival area will be absorbed by mitochondria cells in alveolar structure, thus creating an increase in cell metabolism, resulting in shortened time for orthodontic treatment.

It might, nonetheless, create a warm sensation/tingling in the affected area, which is common and normal.

## **Pain Relief Caused by Dental or Orthodontic Treatment**

The device supports the healing function of the cells and reduces the pain that might be caused by dental or orthodontic treatment in the jaw or tooth area.

## **Safety and Environmental Considerations**

Intense near-infrared light, emitted from LEDs (Light Emitting Diodes), is beneficial, and no side effects have been reported for using healing light in the near infrared spectrum.

## **Limitation for Use**

In case of the below conditions and contraindications, the use of this device is NOT recommended at all.

- Photosensitivity or use of drugs that may cause photosensitivity
- Poor oral hygiene
- Use of osteoporosis drugs
- Epilepsy
- Periodontal disease or acute oral infection
- Patients who has implanted cardiac device unless that device is known as it would not be affected by magnetic fields

## PRECAUTIONS

### Usage

- Make sure that you are entirely familiar with the correct procedures for operating the device prior to the usage.
- Stop using the PBM Ortho device or its accessories if you have an allergic reaction to it. If this happens, contact your doctor or orthodontist immediately.
- Do not stare directly & closely at the light source on the mouthpiece; it might cause eye irritation especially for children.
- Do not use the device near high frequency (HF) equipment.
- Although it can be used while lying down, it is recommended to use while sitting.
- There is no need to be motionless during the use; you might watch TV, read, etc.
- Do not use the device while performing complex tasks or operating machinery.

### Charging

- During charging, it is not recommended to do a treatment session or have the mouthpiece attached to the battery pack. When the charging period is completed, then you can use the device freely.
- During charging, please position the battery on a stable flat surface and out of the way to avoid tripping hazards.
- The battery is not water resistant and should not be rinsed or submerged in water.
- Do not incinerate the battery by exposing it to excessive heat, short circuit or any similar hazardous action to the battery.
- Mishandling the battery may result in a burn, fire, or possibly even an explosion.

### Mouthpiece

- During usage, you should bite gently on the bite pad of the mouthpiece. Chewing or clenching on the mouthpiece may harm the device, or lead to a choking hazard.
- The mouthpiece must be detached from the battery when cleaning and should be kept in a dry environment.
- The mouthpiece is not dishwasher safe.
- The mouthpiece should not be dried with a hairdryer
- Exposing the mouthpiece to excessive heat and elevated temperatures may result in the softening of its material, which would be considered permanent damage of the device.

### **Additional Precautions**

Do not attempt to repair the device or battery. No servicing or modification of this device is permitted.

PBM Ortho has been tested and met the requirements to provide reasonable protection against Radio Frequency Interference in a residential setting. However, there is no guarantee that interference will not happen in a particular situation.

#### **If interference with a radio or television is observed it is recommended to:**

- Connect the equipment to a different outlet.
- Increase the distance between the equipment and the receiver.

### **Additional Precautions for Children**

- Keep the PBM Ortho device away from pets and children.
- It is not a toy and should not be played with.
- When it needs to be used by a child, the steps of usage & sequencing must be clearly demonstrated for them.
- Subsequently, each session of use should be supervised by a parent or an adult to ensure correct usage of the device.

### **Disposal to Environment**

- The User Manual and packaging are recyclable and should be disposed of with other recyclable wastes.
- In order to preserve the environment, the device should not be disposed of with normal household wastes. It must be delivered to the collection points designated for the recycling of electrical and electronic equipment wastes.

## **OPERATING CONDITIONS**

Ambient temperature range: 5°C to 45° C

Maximum temperature on any external component at maximum ambient temperature: 40°C

Relative humidity range: 15 to 93 % non-condensing

Atmospheric pressure range: 700 to 1060 hPa

## **TRANSPORT AND STORAGE CONDITIONS**

Minimum ambient temperature: -20°c

Maximum ambient temperature: 55 ° C

Maximum humidity: 93% non-condensing

## **WARRANTY MANUFACTURER**

The manufacturer will only warrant the original purchases from PBM Healing for flaws in workmanship & material. There is a two (2) year warranty from the date of the original purchase.

## **HELP**

Please contact your doctor or PBM Healing rep for all inquiries regarding your device usage.

Please contact our Support Center for help, if needed, for setting up, using or the device. To report unexpected events or operation and troubleshooting.

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